


EZ9™ and EZ10™

Closed Door System Fully Automatic Autoclaves

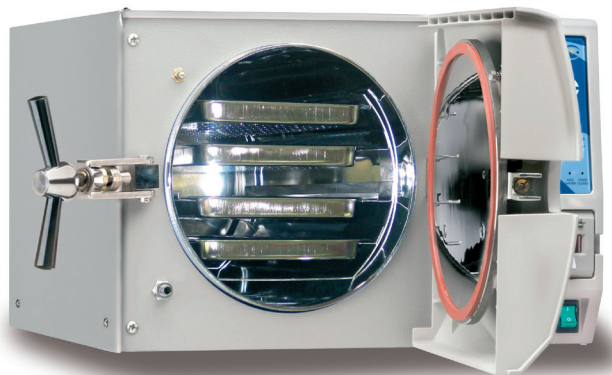

EZ10 with 4 large trays

EZ9™ and EZ10™

Autoclaves will meet all your sterilization needs

- Chamber door remains closed during entire cycle
- Hepa filtered air is circulated in chamber during drying
- All program parameters can be customized and stored
- New brighter and easier to read panel
- Two year parts and labor warranty

* Above photo shown with optional printer

EZ9™ and EZ10™ Closed Door System Fully Automatic Autoclaves

With the simplicity of one touch design all your sterilization and drying needs are fulfilled. The Tuttnauer fully automatic autoclaves fill, sterilize, exhaust and dry at the touch of a button.

EZ9 and EZ10 Series:

Item No.	Chamber Dimensions Ø x L In(mm)	Volume	Overall Dimensions DxWxH In(mm)	Standard Cassettes	Tray Dimensions DxWxH In(mm)	No. of Trays	Standard Unwrapped Cycle time	Voltage(V) Freq.(Hz)	Power(W) Current(A)	Shipping Weight lbs(kg)
EZ9™	9" x 18" (230 x 470)	5gal/19L	21.5"x20"x14.4" (545x510x365)	2 full/2 half	16.3"x6.7"x0.8" (415x170x20)	3	Cold: 27 mins Hot: 13 mins	120V 50/60	1400W 12A	95 (43)
EZ10™	10" x 19" (254 x 475)	6gal/23L	21.5"x20"x14.4" (545x510x365)	3 full/3 half	16.3"x6.7"x0.8" (415x170x20)	4	Cold: 30 mins Hot: 14 mins	120V 50/60	1400W 12A	95 (43)
*EZ10k™ Kwiklave™	10" x 19" (254 x 475)	6gal/23L	21.5"x20"x14.4" (545x510x365)	3 full/3 half	16.3"x6.7"x0.8" (415x170x20)	4	Cold: 21 mins Hot: 11 mins	230V 50/60	2200W 10A	95 (43)

*Model EZ10k™ requires a voltage between 220V and 235V for proper operation. When applicable a buck/boost transformer (0.5KVA) may be required.

Features:

- Closed door active drying system to maintain sterility and ensure efficient drying of packs and pouches
- High efficiency air pump circulates air in chamber during drying cycle
- 0.2µm Hepa air filter (replaceable)
- New brighter and easier to read panel
- Program parameters can be customized and stored

Control Panel:



1. Program Selection
2. Detailed Digital Display
3. Programmable Keys
4. System Indicators


Optional Internal Printer



Documents:
date, time,
temperature,
pressure.



Clean&Simple™ Ultrasonic Cleaners and Tablets

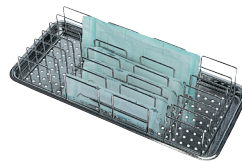

Instrument Pouch Rack



Steam Distiller



Chamber Brite Autoclave Cleaner



Tuttnauer USA Co., Ltd.

25 Power Drive, Hauppauge, NY 11788 Tel: +800 624 5836, +631 737 4850 Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com, www.tuttnauerUSA.com

Tuttnauer
Your Sterilization & Infection Control Partners